

# CAREpoint EMS Workstation

## Application Note: 183V020A

**Title:** How to configure Windows Update Settings

### Overview:


This procedure guides you through the steps for configuring Windows Update Settings. This procedure will also guide you through disabling an update that Windows will continue to download and install over and over again. This procedure should be performed by someone who is familiar with the CAREpoint Workstation.

### Execution:


1. In the lower right corner of the screen, double click on the yellow Windows Update shield.



2. An Automatic Updates window will appear, pick “Custom Install (Advanced)”, and then click Next.

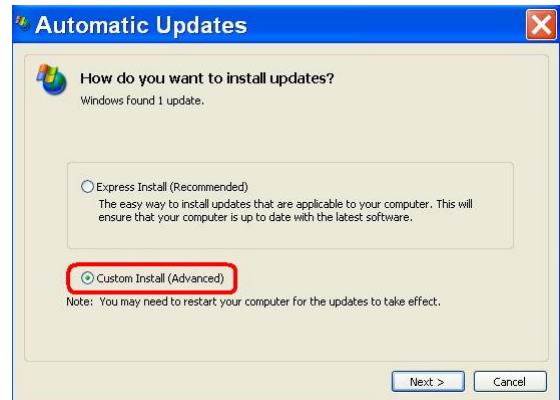

3. Locate and uncheck the checkbox that says “Security Update for Windows XP (KB938464)”. Leave any other checkboxes checked. Then click Close.

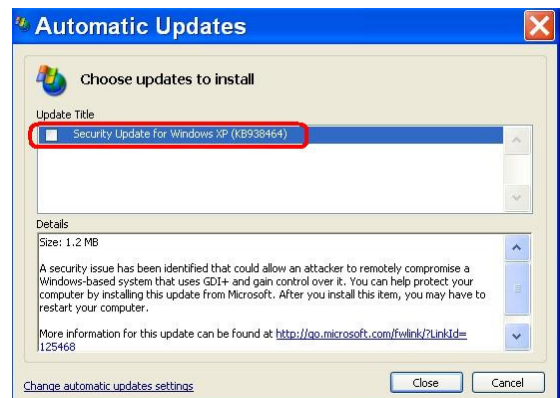


4. A Hide Updates window will appear. Check the checkbox that says “Don’t notify me about these updates again.” Then click OK.



5. On the blue, on-screen keypad, slowly type #-E-X-P-L-O-R-E (#-3-9-7-5-6-7-3). This will open a Windows Explorer window.



6. On the left side of the Windows Explorer window, right click on “My Computer” and then click on “Properties”.

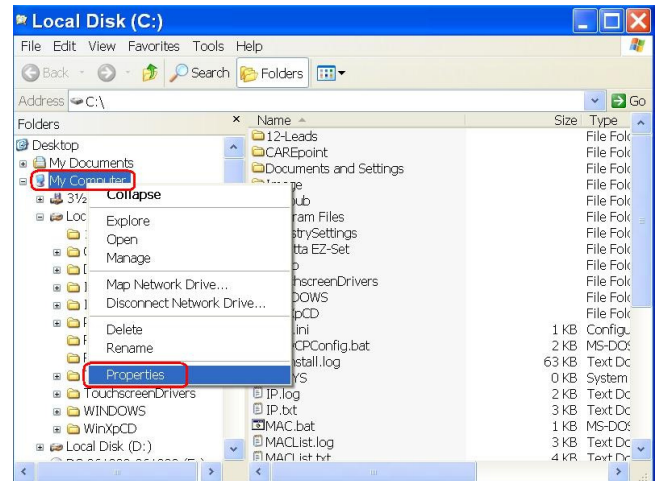


7. A System Properties window will appear. Click the “Automatic Updates” tab at the top.



(Continued on next page)

- Pick “Download updates for me, but let me choose when to install them.” Then click OK. Click the “X” in the upper right corner of the Local Disk (C:) window. This will bring the main CAREpoint screen back up.

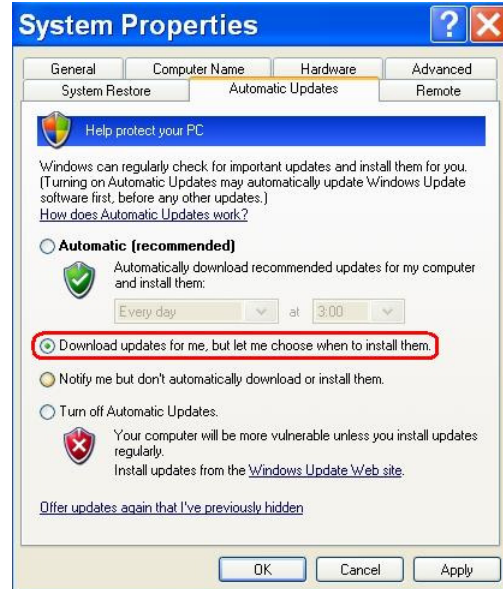

**- END PROCEDURE -**